



Instruction for 'BRIS' Scanner software installation on 'Windows10' computer.


1. Insert the memory stick that comes in a kit with the scanner to the USB connector.

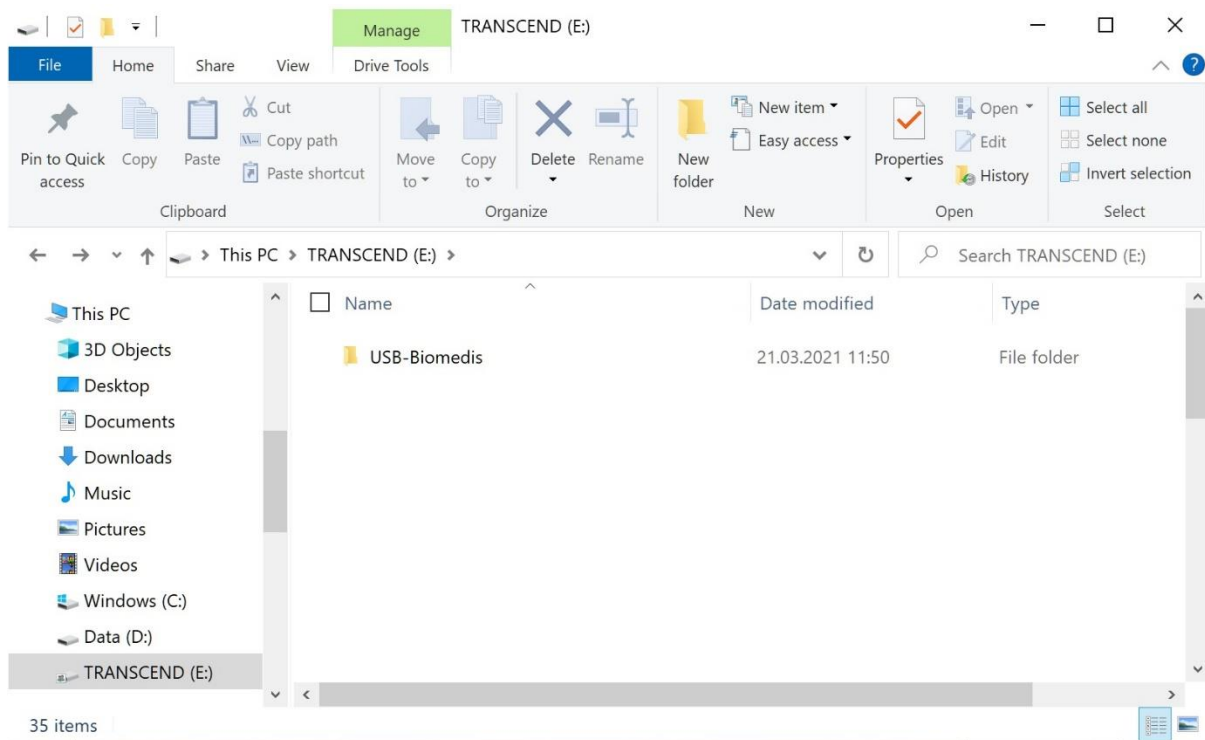


2. Open 'USB-Biomedis' folder:

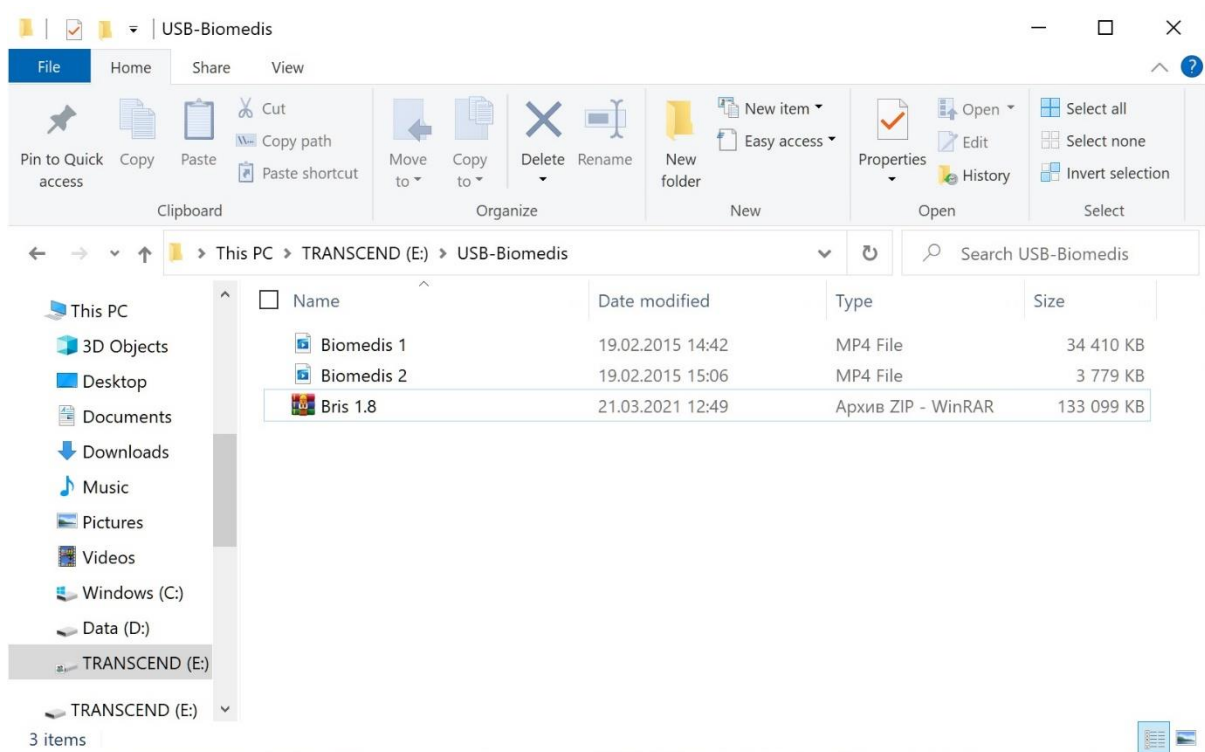


3. Create on disk C: folder 2002

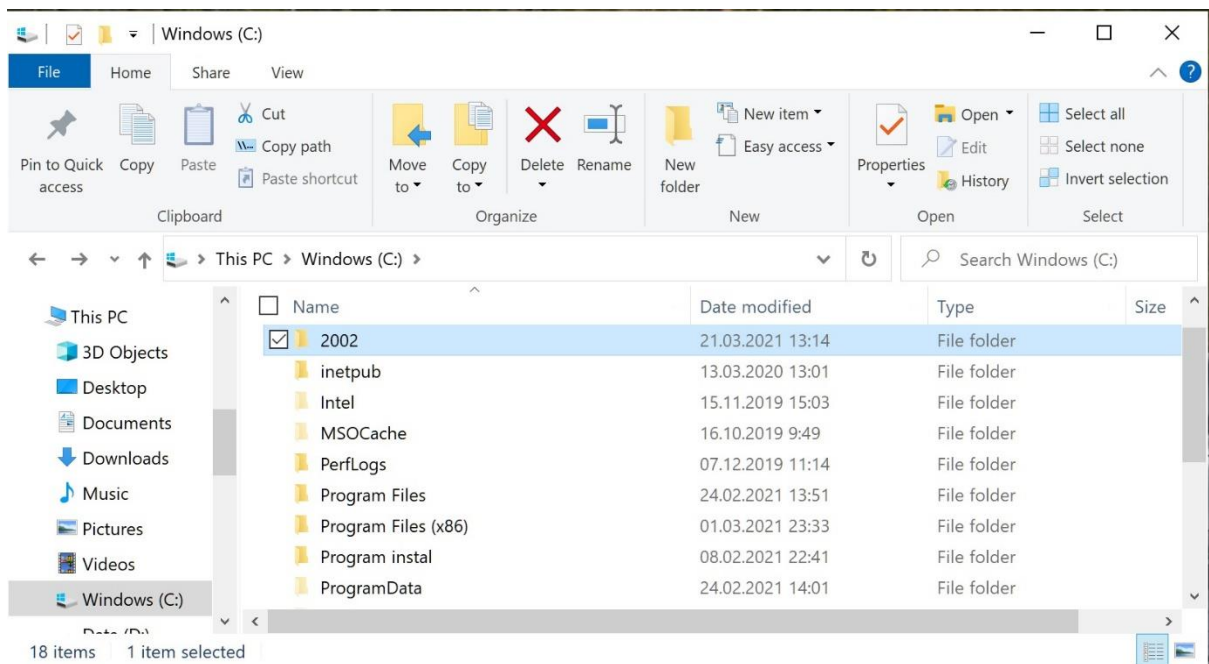


4. Unpack the archive "C:\2002\" into the folder 2002

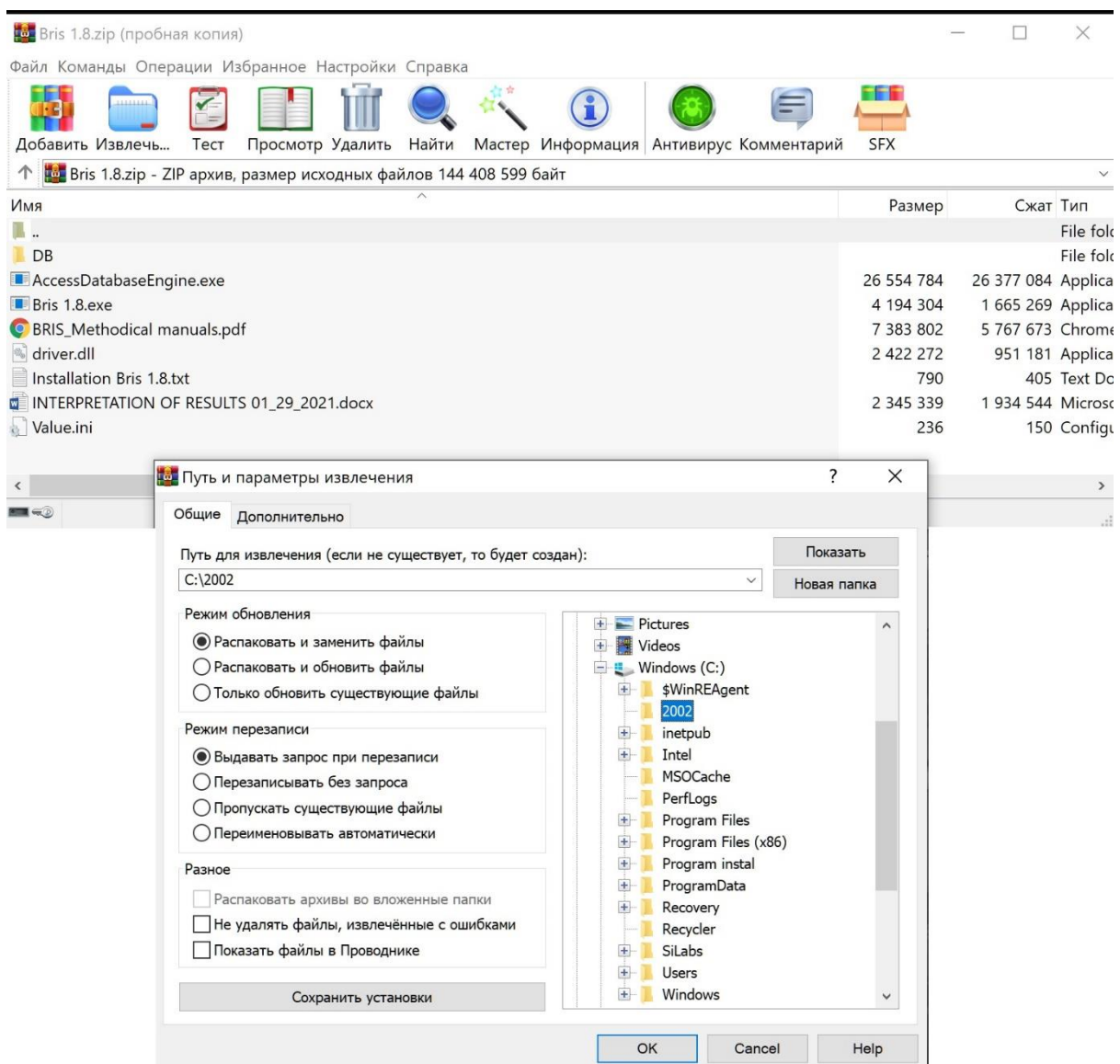

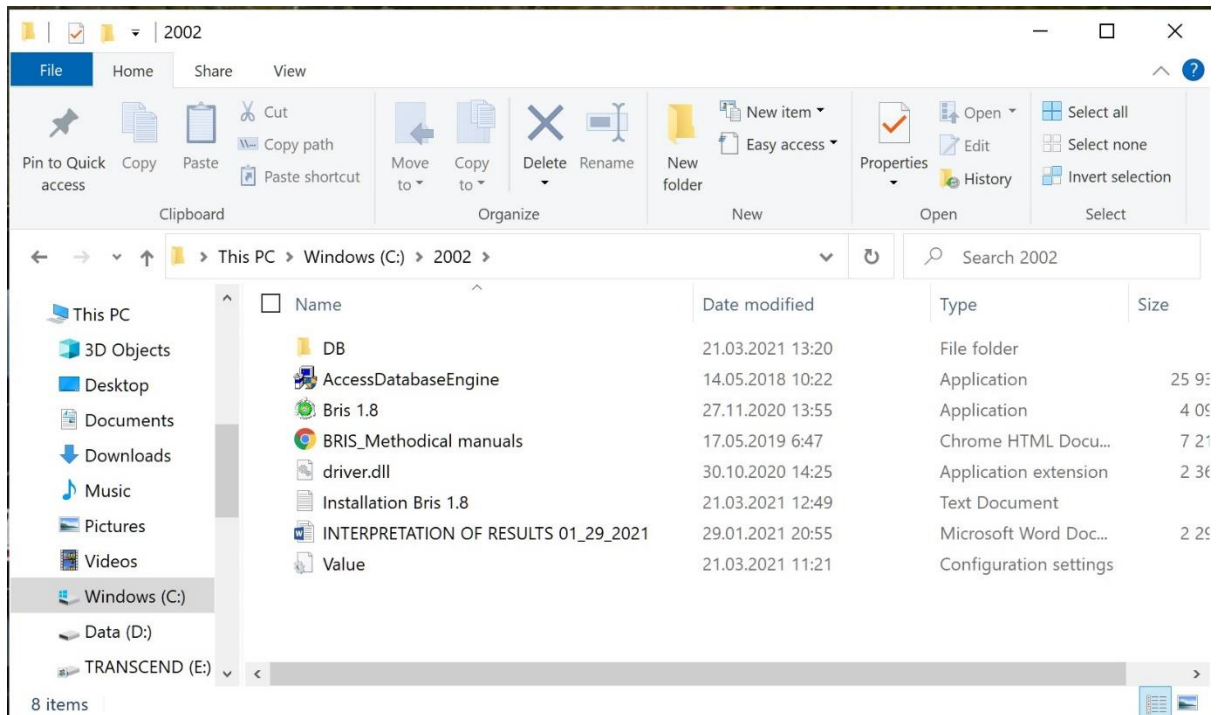


5. Put the tag mark of the 'Bris 1.8' program to the desktop.

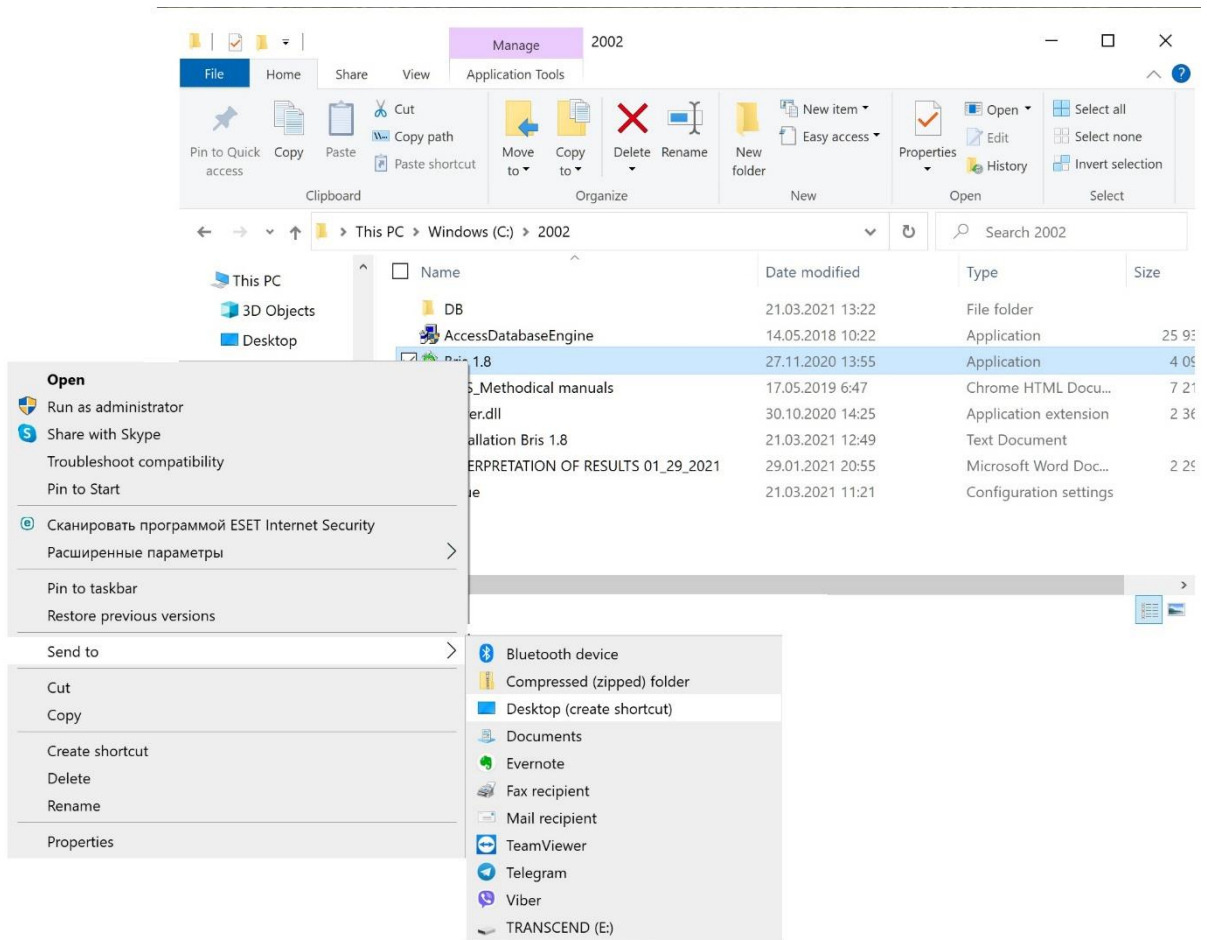


6. Launch the program “Bris” by pressing the icon




When launching the program will be started in English. For correct registration of the users in the program it is needed to complete the step by step setting, described in points: 7,8,9,10.

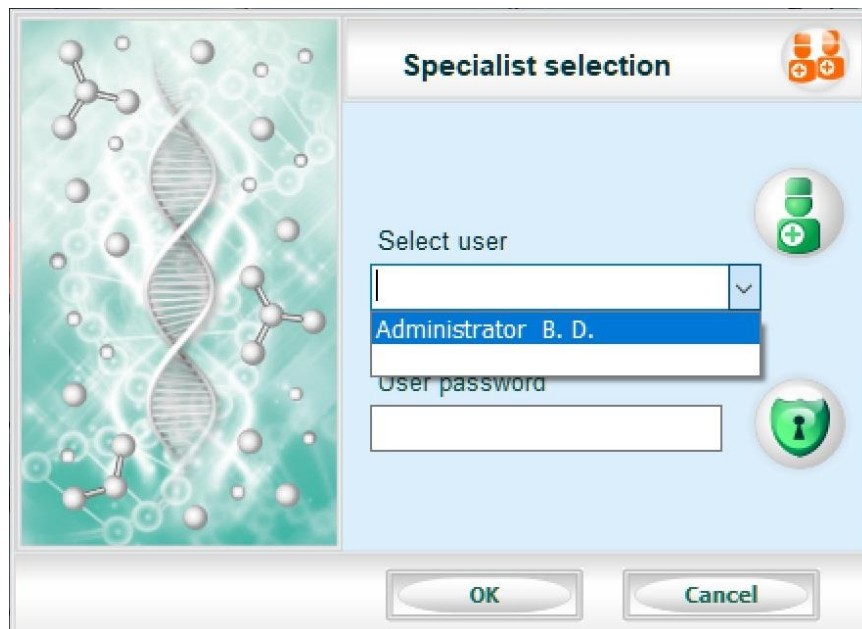


7. To change the tab page ‘Format’ into Russian it is needed to open ‘Control panel’ and select tab page ‘Time, Language, Region’ or ‘Regional standards’

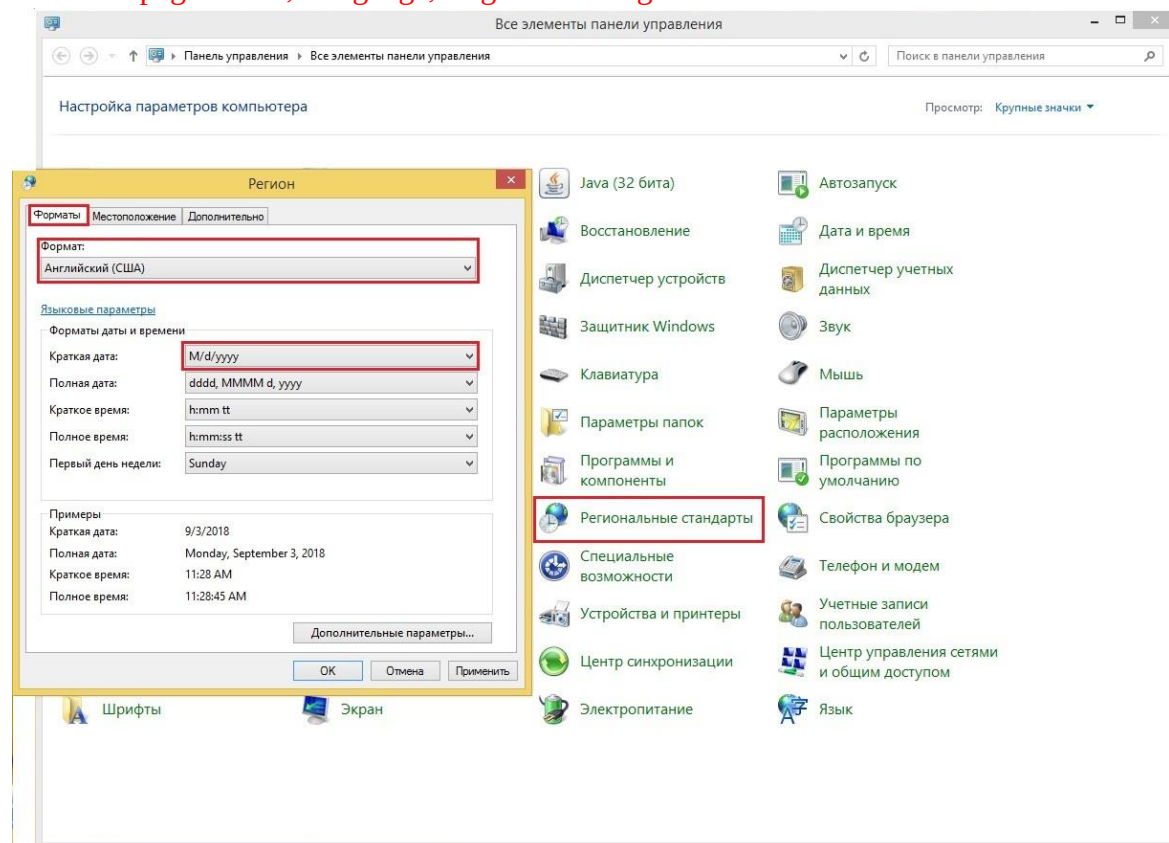


8. Select the ‘Russian region’ format in the tab page ‘Format’ and install it.

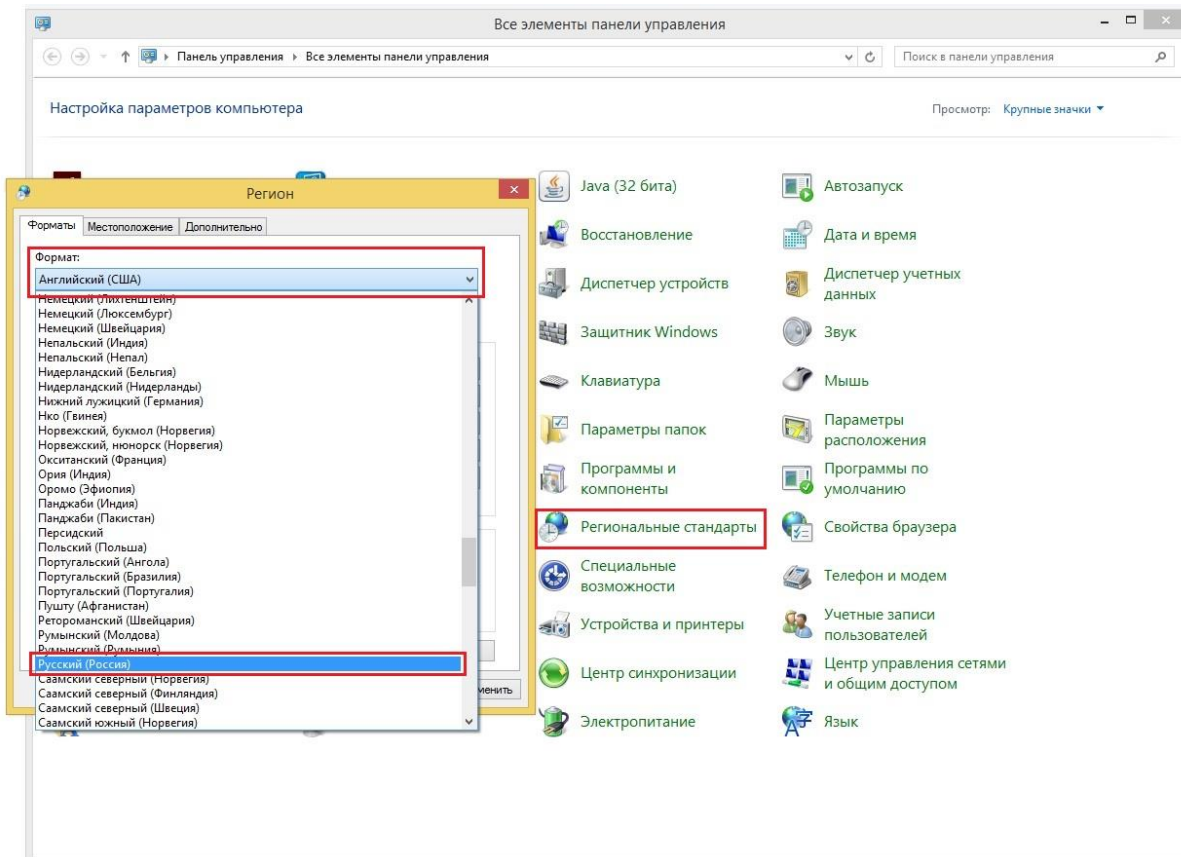

9. In the window of 'Format of date and time' in the section 'Short date' instead of (dd/MM/yy) select (dd.MM.yyyy):

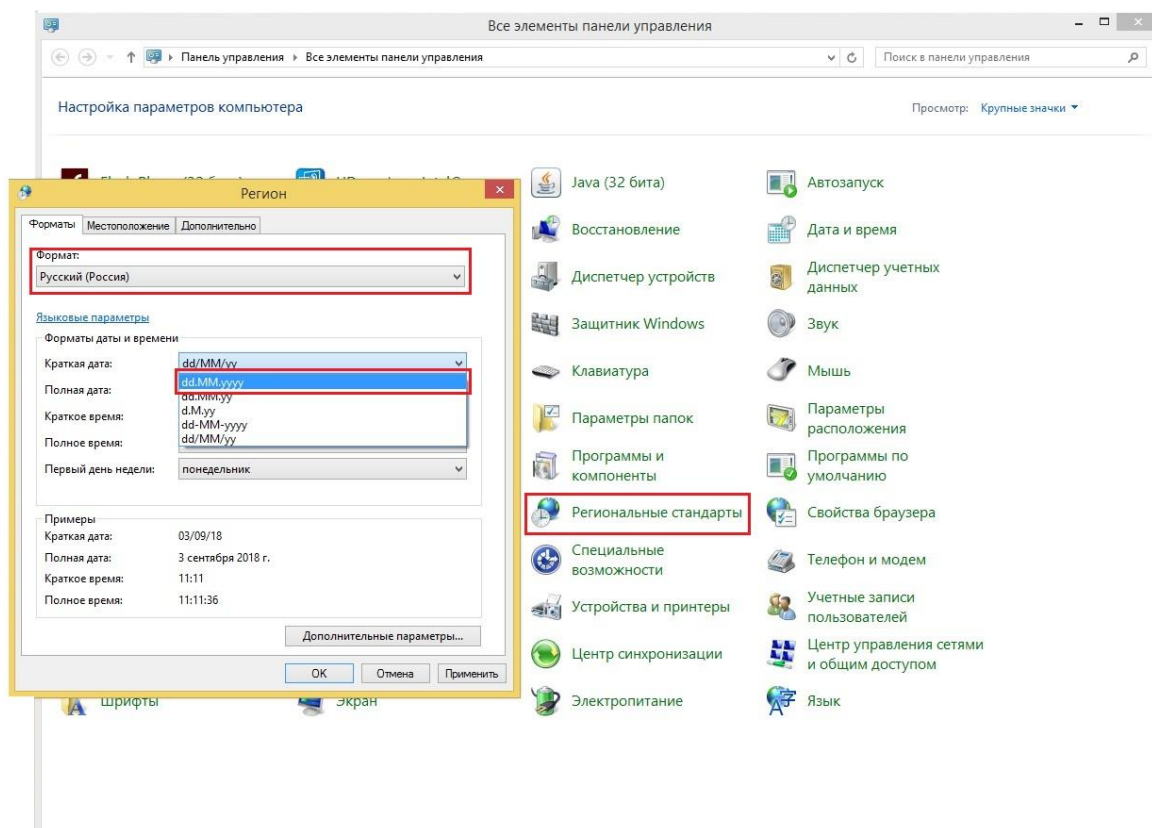


10. Нажать «Применить» и далее «ОК»

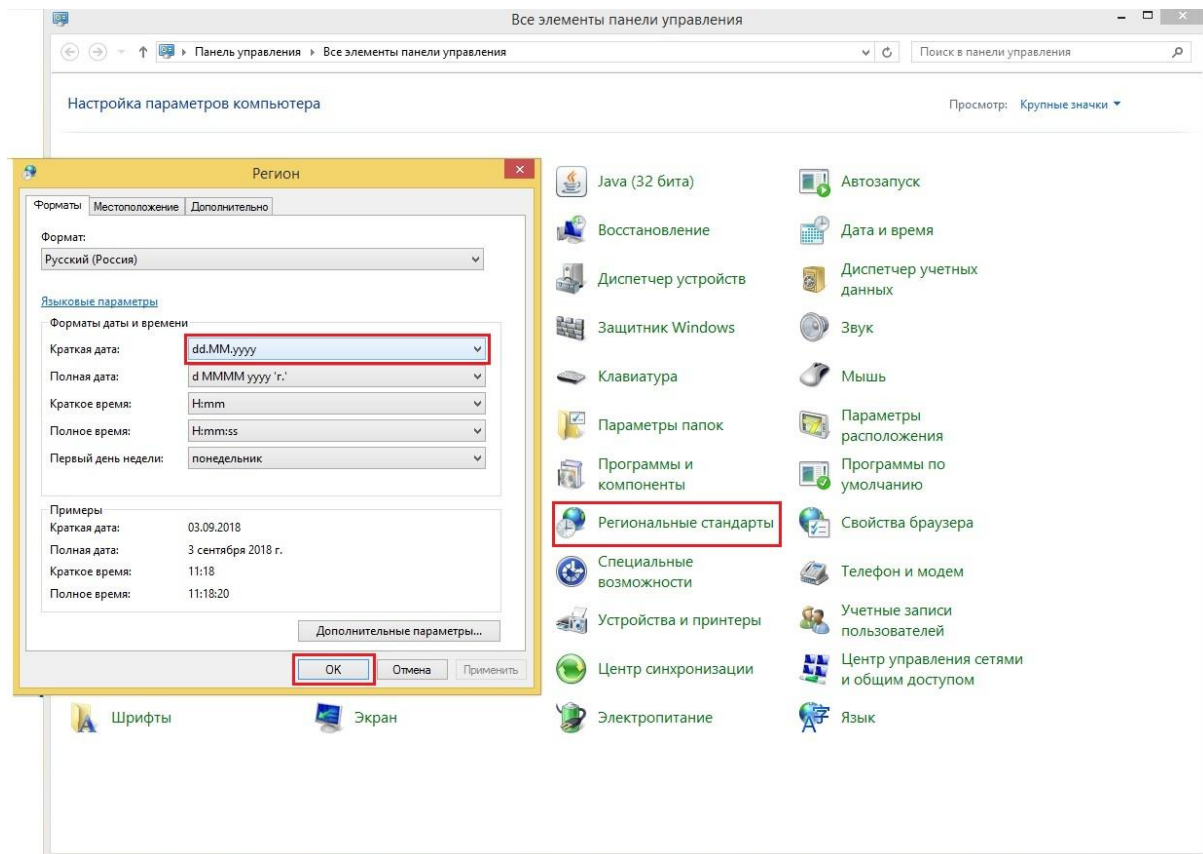


11. Relaunch the program. It will be started in English.

12. Select «Setup».

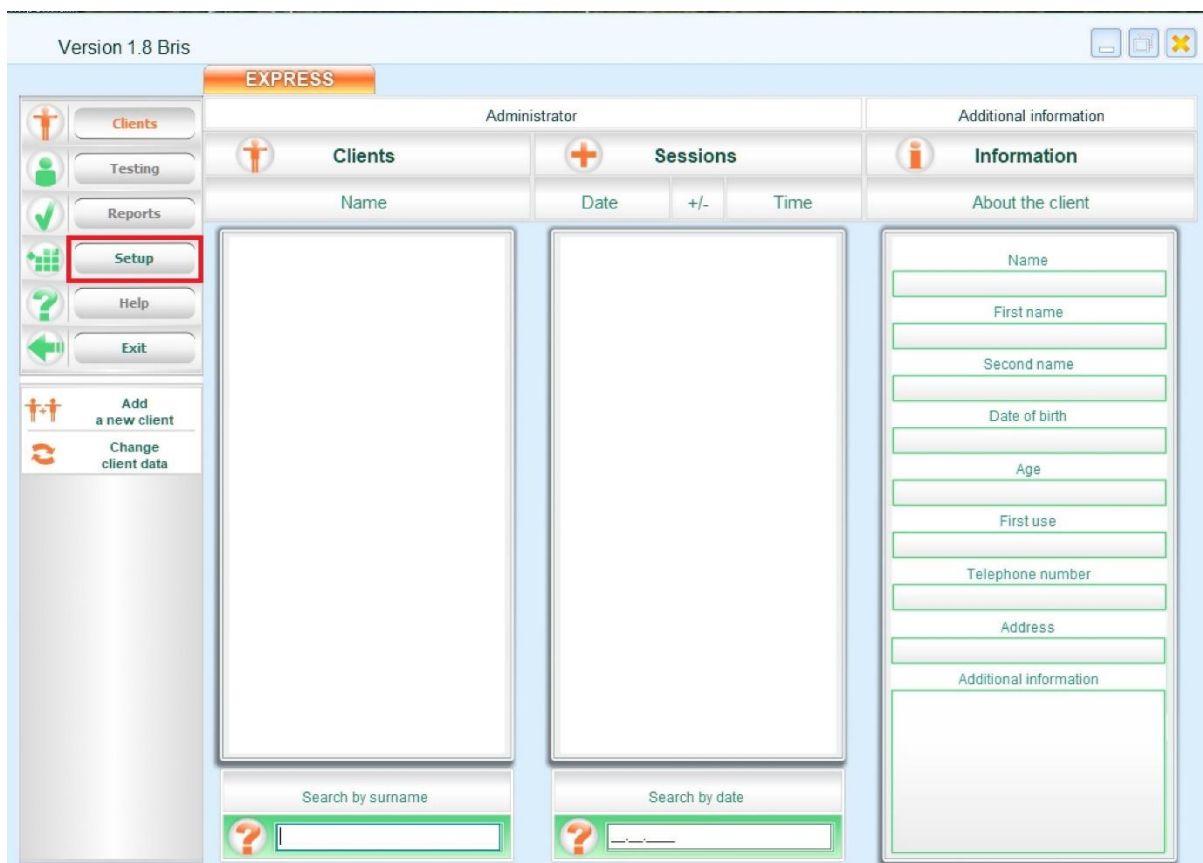


13. Select the tab page 'Device Registration' and in the opened window text the registration code of the device. The registration code can be found on the back side of the device or in the passport of the device.


ATTENTION, EACH DEVICE HAS ITS INDIVIDUAL CODE:

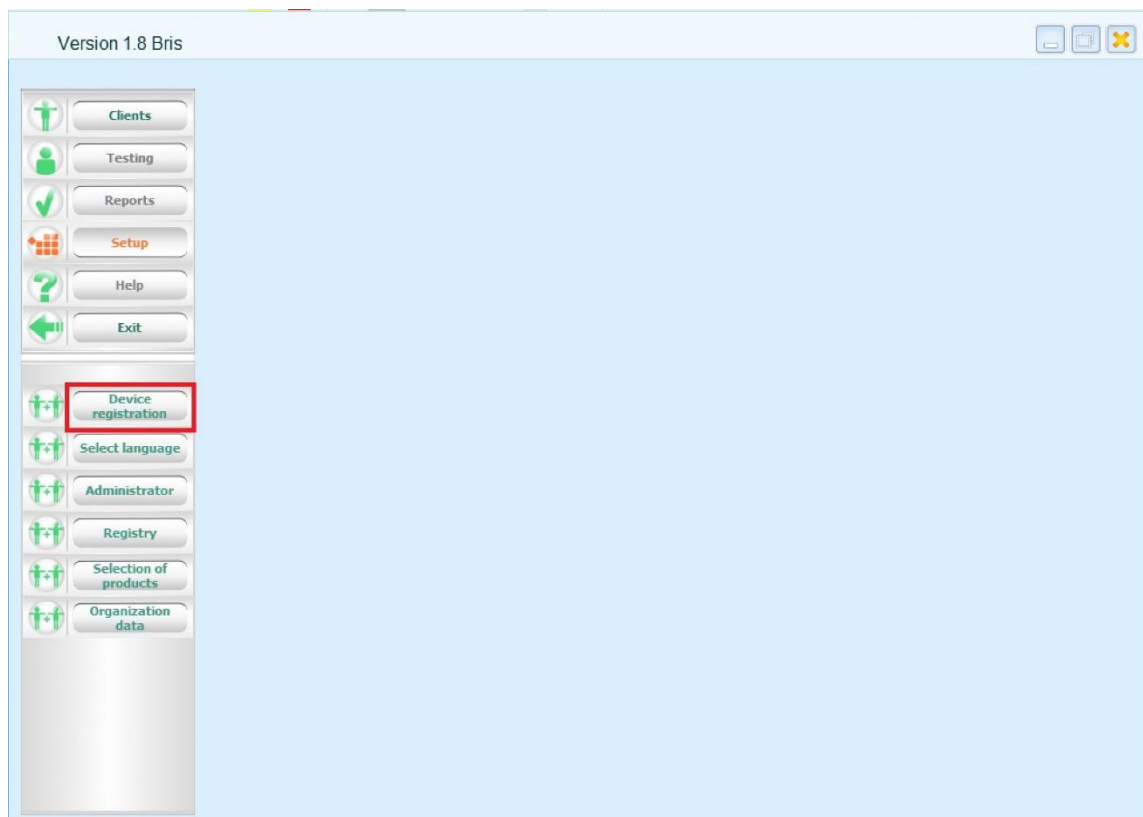


(example of the device registration is in the following picture)

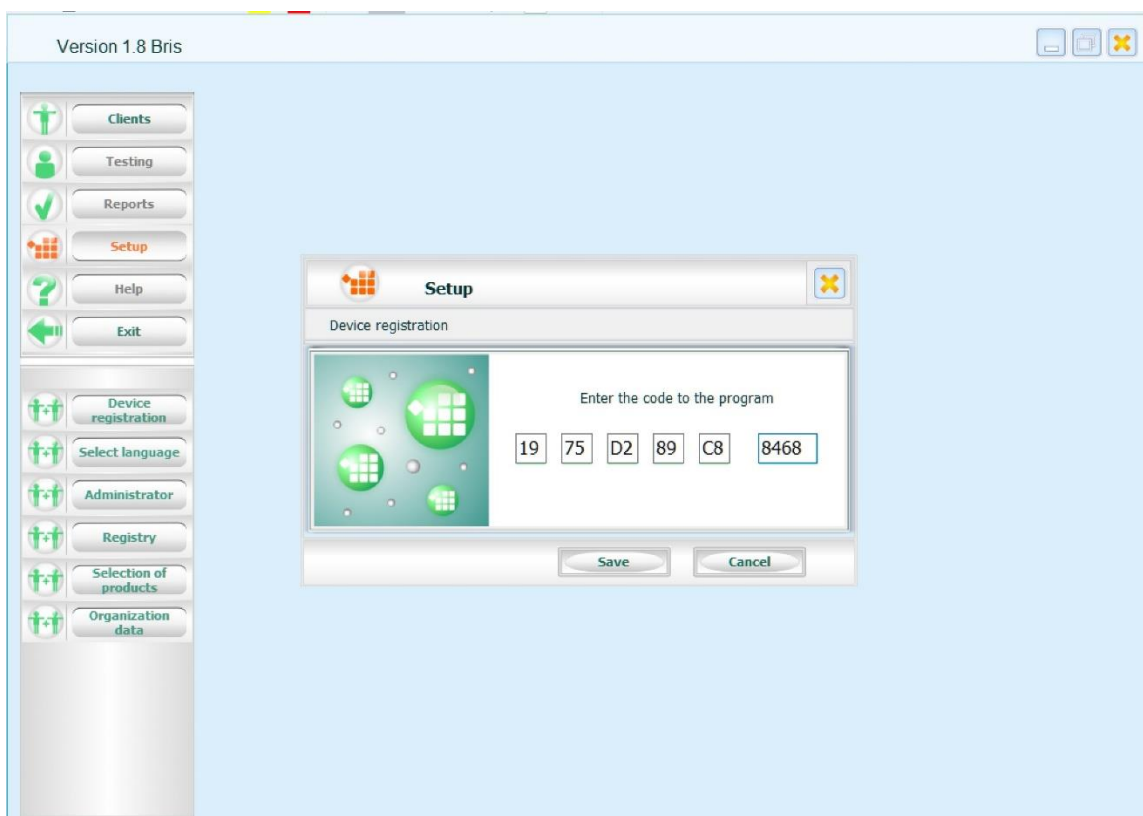


14. As soon as you enter the serial number of the device, press 'Save' and the program will remember your device.

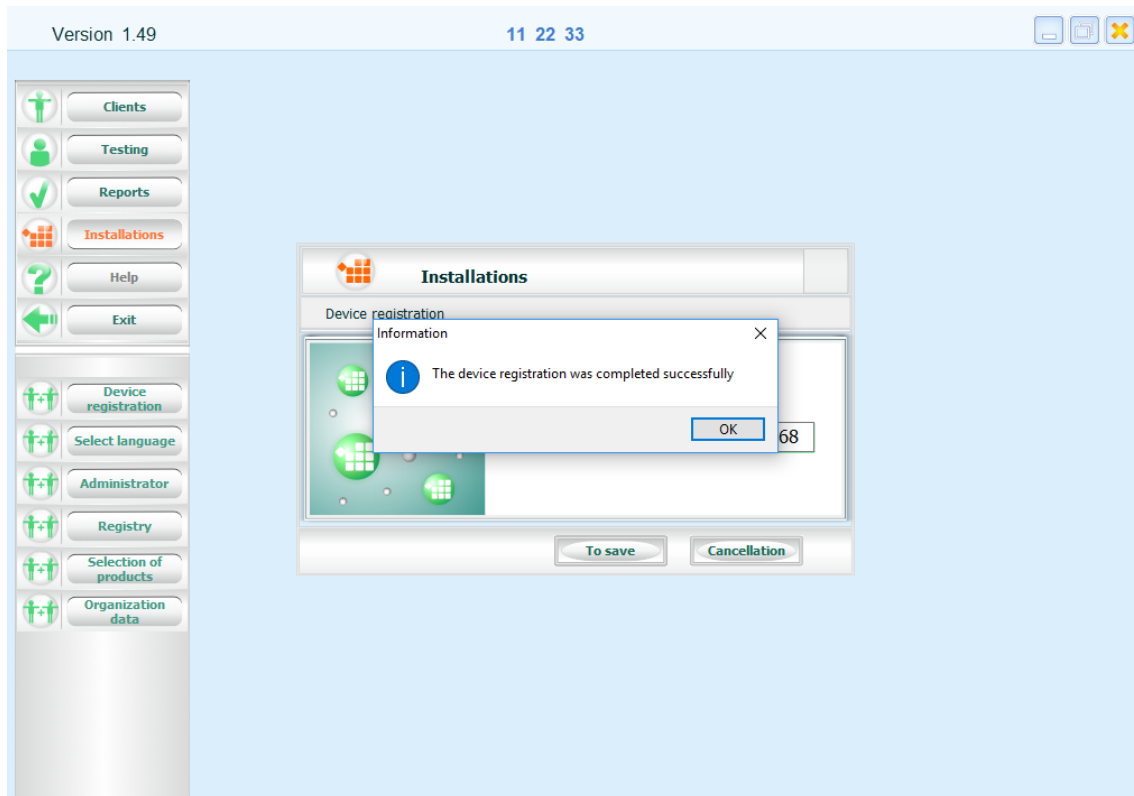


15. As soon as you made all the steps you can change the language in the 'Installations' by selecting the tab page 'Select Language'. In the opened window you select the interested language and press 'To save' (currently there are no translations to Turkish and Arabian languages).

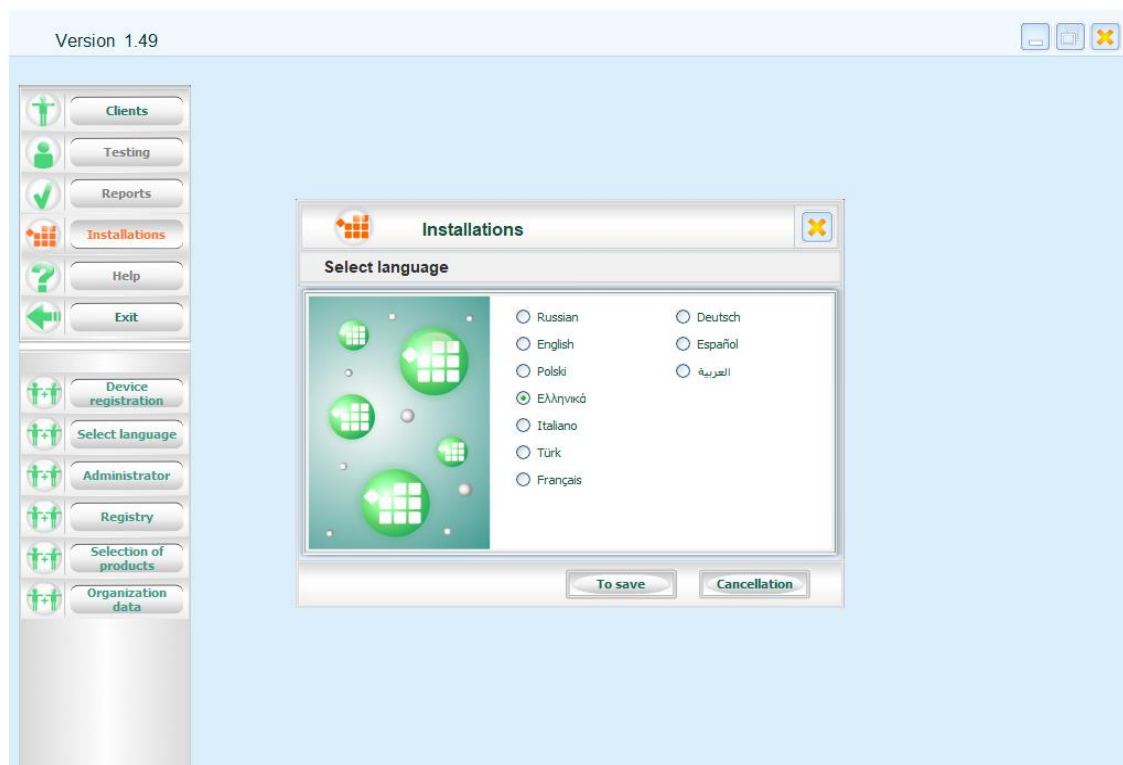


16. After selecting the language the program will request the restart. Exit the program and launch it again.

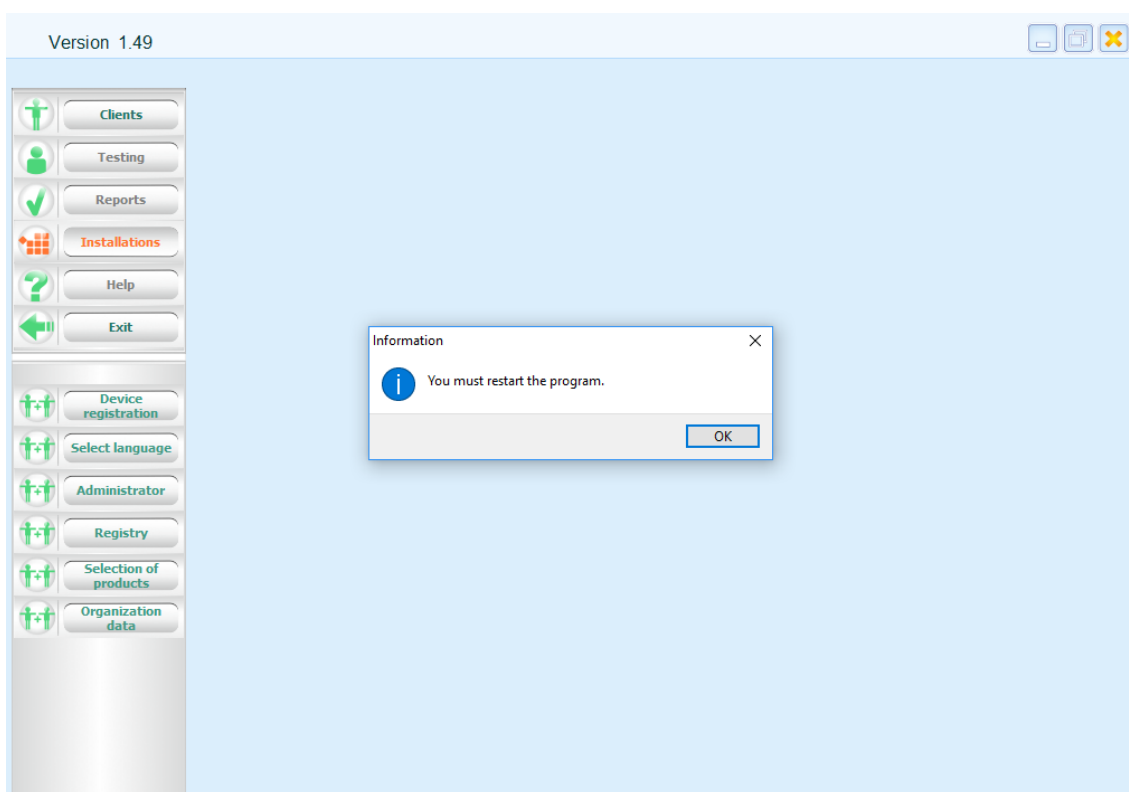


Then the program will be launched in the selected language (the example is in Greek)

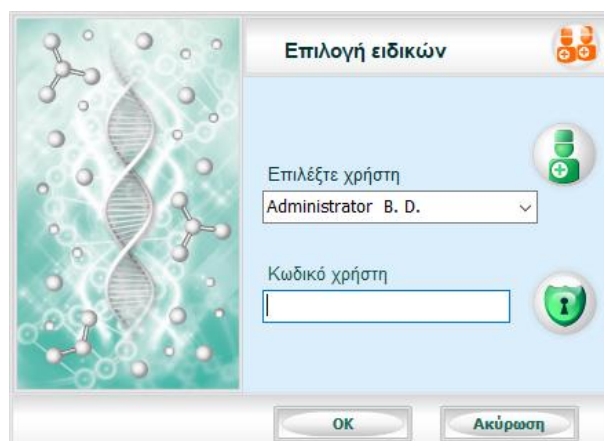

ON COMPUTERS when trying to enter the program there can appear the following message:

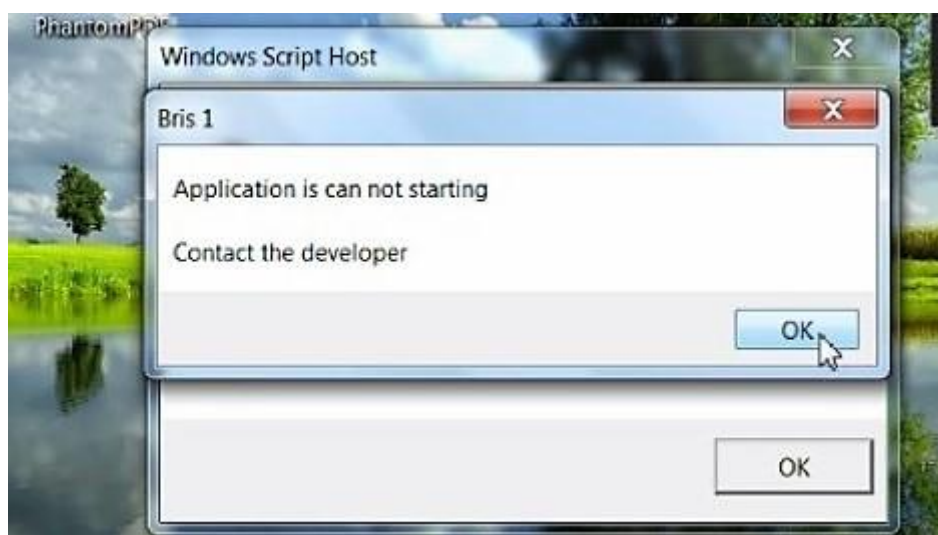

To solve this problem it is needed to open disk C, open the '2002' folder and launch the installation of the 'AccessDatabaseEngine' program.

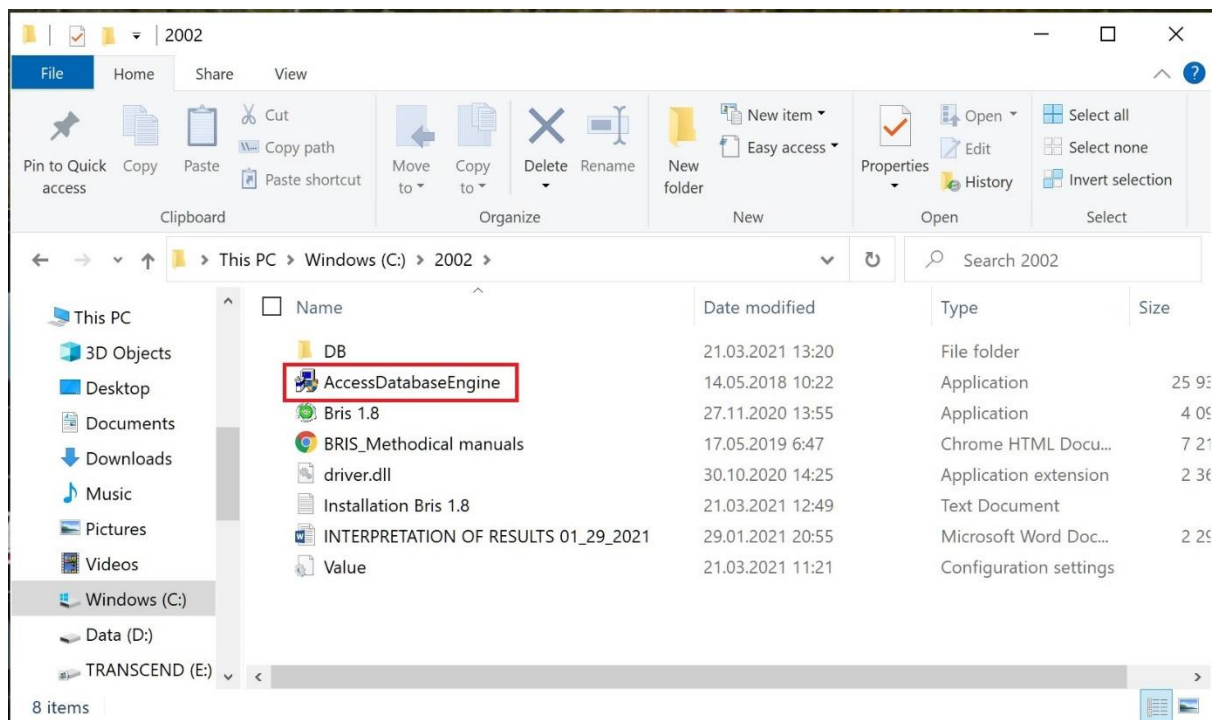

After these steps launch the icon 'BRIS' on the desktop of your computer.



Use the program for the benefit of yourself and others!